

## Charisma® EasyShine

### High lustre in one single step.

Has polishing ever been so easy?

- » **Excellent results:** The diamond-enhanced formula delivers excellent long-lasting polishing results through a pressure controlled polishing process.
- » **One-step polishing system:** No need to change instruments. Smoothing and high-lustre polishing in one step are possible by adapting the applied pressure. That saves your time.
- » **Economic solution:** It is ideally suited for all Charisma composites, such as Charisma Diamond, Charisma Topaz, Charisma Classic. All polishers are suitable for sterilisation and therefore reusable.

Giving a hand to oral health.



**KULZER**  
MITSUI CHEMICALS GROUP



Filename: HER1706_11534_Charisma_Easy_Shine_Sellsheet_EN_Kulzer	Kulzer Toolbox No: W18834	Size closed (EU, USA): EU 210x297mm	Number of pages 2 pages
SAP No.: -	ORT-Version: 02	Speciality in printing/Spotcolours: -	Lacquer: no
Marketing Manager: Sebastian Oyegbesan	Producer: ORT	printing profile: iso coated v2 (eci)	Date 08-03-2017 Kulzer final release 08-17-2017

# Charisma® EasyShine

High lustre in one single step.

Equipped with the well-known Charisma quality for excellent polishing results:

- **Fast and easy:** The diamond-enhanced formula delivers excellent polishing results through a pressure-controlled polishing process. Highly polished surfaces prevent plaque accumulation and reduce the risk of secondary caries for long-lasting natural restorations.
- **One single step:** Smoothing and high-lustre polishing in one step are possible by adapting the applied pressure. There is no need to change instruments which saves time and brings less technique sensitivity.
- **Universally applicable:** The three available shapes (Flame, Large Flame and Cup) cover all intra-oral applications. They are compatible with all kinds of composites and are ideally suited for all Charisma composites from Kulzer.
- **Durable and reusable:** All polishers are suitable for sterilisation and therefore reusable.

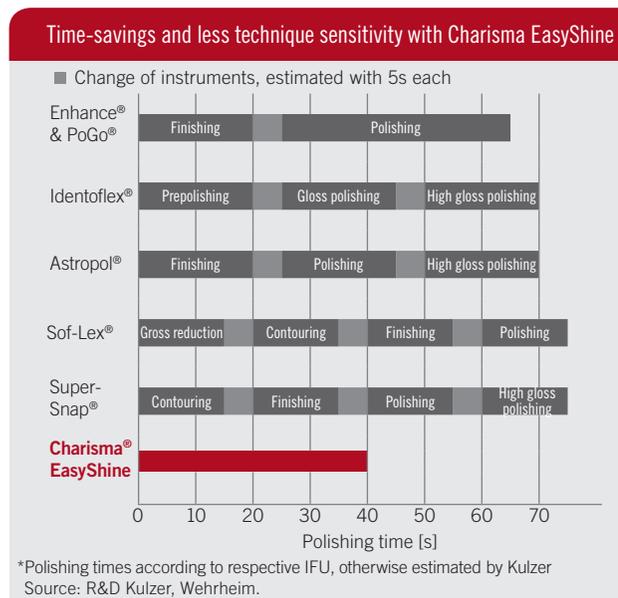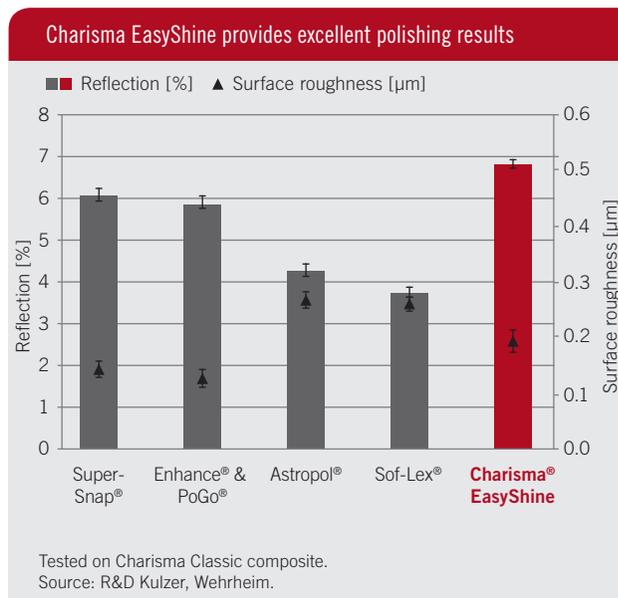

Charisma EasyShine is made of a polyurethane matrix which contains an abrasive mixture of microfine diamond powder and silicone carbide. It is that premium combination of high-quality abrasive and specifically modified bonding which allows smoothing and high-lustre polishing in one step, just by adapting the applied pressure. This pressure-controlled polishing process saves your time and enhances your restoration results.

Smoothing and high-lustre polishing are done without the need to change instruments. Polishing in one single step reduces the polishing time and is less technique sensitive.

Charisma EasyShine makes your restorations shine. The benefits are yours to max out.

Art. code	Article	Content
66060002	Charisma® EasyShine Polishing Kit	2 x shape Flame, 2 x shape Large Flame, 2 x shape Cup
66060003	Charisma® EasyShine Refill – Flame	12 x shape Flame
66060004	Charisma® EasyShine Refill – Large Flame	12 x shape Large Flame
66060005	Charisma® EasyShine Refill – Cup	12 x shape Cup



Super-Snap® is a trademark of Shofu Dental GmbH · Enhance® and PoGo® are trademarks of Dentsply De Tray GmbH · Astropol® is a trademark of Ivoclar Vivadent AG · Sof-Lex® is a trademark of 3M Company · Identoflex® is a trademark of Kerr Corporation

© 2017 Kulzer GmbH. All Rights Reserved.

Kulzer GmbH · Leipziger Straße 2 · 63450 Hanau · Germany · kulzer.com

W18834 EN 07/2017 ORT