


Ivory® Rubber Dam Clamps

A secure hold for absolute isolation.

Choose from a wide variety of uniquely designed Ivory Rubber Dam Clamps to find the perfect one for each tooth. Our line includes a clamp for every application. Every single clamp is crafted in the USA from high-grade stainless steel and carefully hand-set so that you can confidently affix the rubber dam whenever it's needed.

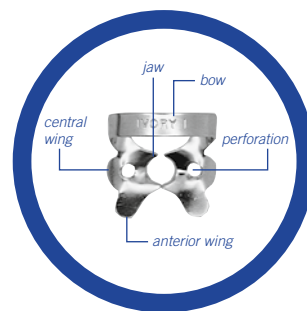

Molar

3



For small lower molars. Narrow, flat jaws so as not to impinge in tissue.

3
W3

4



For small uppers. Narrow jaws adapt to fit round cervical (neck) collars.

4
W4

5



For large upper molars. Jaws adapt to curvature at the gum line.

5
W5

7



General purpose lower molar. Medium-sized bow and jaws. Flat jaws. For teeth which tend to be square in form.

7
W7

7A



For large lower molars. Same as No. 7, only larger. Large bow, broad jaws.

7A

8



For small upper right and lower left molars with short crowns. Buccal jaw broader than lingual.

8
W8

8A



For upper, lower molar roots. Four-point jaws grasp root at corners. Jaws set dam beneath gum line when released.

8A
W8A

10



For small upper left and lower right molars with short crowns. Buccal jaw broader than lingual.

10

11



For small upper right and lower left molars with short crowns. Buccal jaw broader than lingual.

11

12A



For all lower right molars, particularly third molars and partially erupted teeth. Available with serrated jaws only.

12A

13A



For all lower left molars, particularly third molars and partially erupted teeth. Available with serrated jaws only.

13A

14



For upper and lower molars, partially erupted, round teeth or for teeth sloping from the gum line. Small bow.

14
14DT
14T
W14

14A



Same as No. 14, only larger in every way.

14A
W14A

26



For small lower molars.

26

26N



Universal upper with festooned jaws.

26N

27



For large upper and lower molar distal cavities.

27

28N



General purpose lower molar clamp.

28N

56



For all large upper and lower molars.

56
56T
W56

Central Anteriors

00



For small lower anteriors and teeth set irregularly in arch. High bow, narrow jaws.

00

6



For anterior labial and premolar teeth.

6


9



For anterior labial teeth.

9
9T

Tooth Diagram

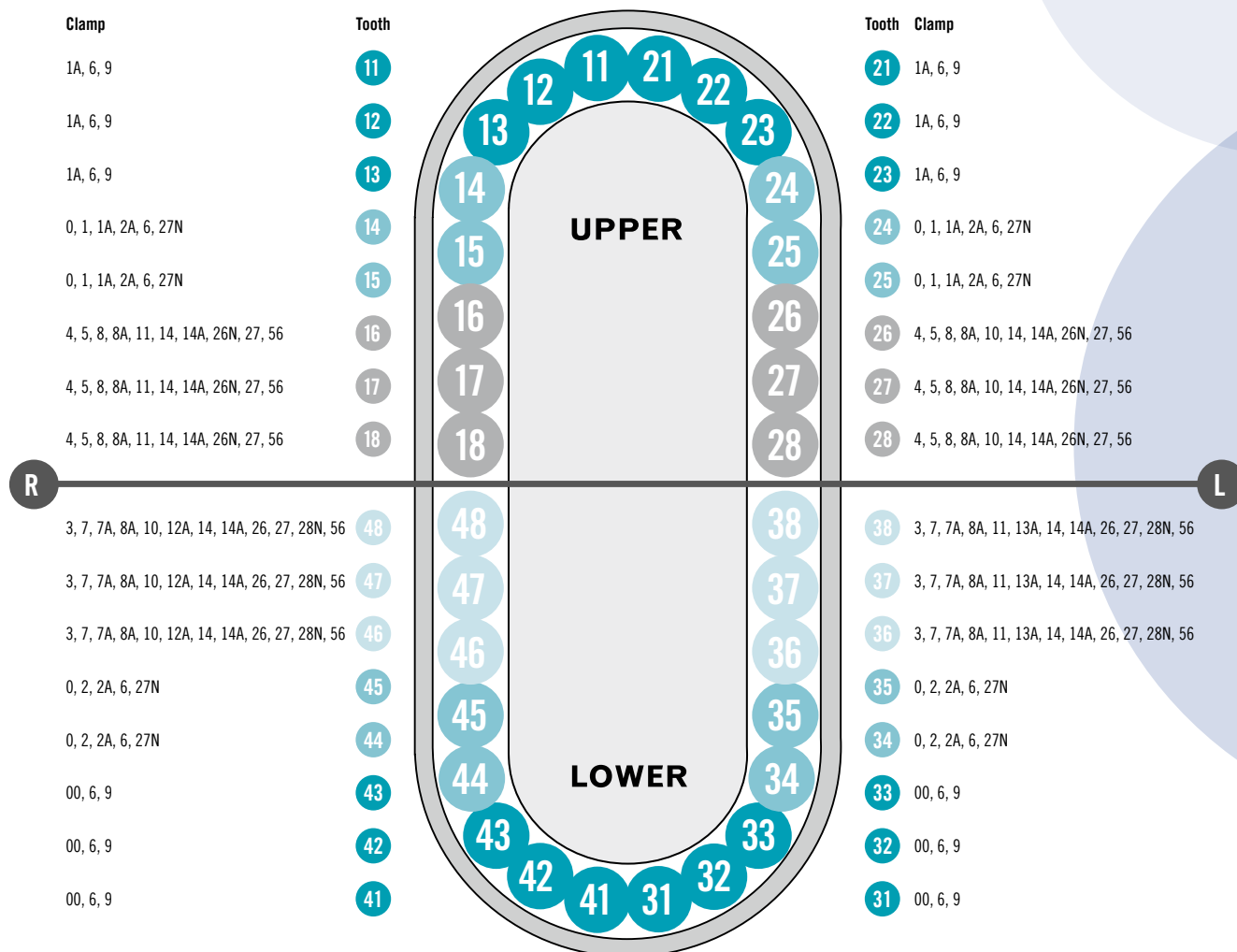

Premolar

0 For upper and lower premolar teeth that are elongated or have small necks. 0 W0	1 General purpose upper premolar. Jaws adapt to curvature at the gum line. 1 1T W1	1A Upper premolar and anterior root clamp. Jaws adapt in a downward slope to grasp roots. 1A
2 General purpose lower premolar. Same as No. 1, only jaws are flat so as not to impinge on tissue. 2 2T W2	2A For upper and lower premolars with large necks. Broad jaws. 2A 2AT W2A	27N Universal premolar clamp. Jaws adapt to curvature at the gum line. 27N

DT Distal access, serrated jaws
T Also available with serrated jaws (Tiger clamps)
A Slightly larger
W/N Wingless

Which Ivory® Rubber Dam Clamp do you need?

Listed below is the complete Ivory Rubber Dam Clamp line, along with guidelines to help you select the right clamp for the job.



Most Popular Clamps

PREMOLARS – WINGED

- 00 Small lower, anterior
- 0 Elongated small neck
- 1 General purpose upper
- 2 General purpose lower
- 2A Large neck premolar

PREMOLARS – WINGLESS

- W2 General purpose lower
- W2A Large neck premolar

MOLAR – WINGED

- 3 Small lower
- 4 Small upper for round teeth
- 7 General purpose lower
- 8 General purpose upper
- 8A Upper and lower molar root
- 12A Right partially erupted, 3rd molar
- 13A Left partially erupted, 3rd molar
- 14 Partially erupted, most popular clamp

MOLAR – WINGED (CONTINUED)

- 14A Large, partially erupted
- 56 All large molars

MOLAR – WINGLESS

- W3 Small lowers
- W8A Upper and lower molar root
- W14A Large partially erupted

Specialty Clamps

PREMOLARS – WINGED

- 1A Root clamp

PREMOLARS – WINGLESS

- W0 Elongated small neck
- W1 General purpose upper
- 27N Universal

MOLAR – WINGED

- 5 Large upper
- 7A Large lower
- 10 Right, small short crown

MOLAR – WINGED (CONTINUED)

- 11 Left, small short crown
- 26 Small lower molar
- 27 Distal access, large molar

MOLAR – WINGLESS

- W4 Small upper
- W5 Large upper
- W7 General purpose lower, square teeth
- W8 General purpose upper, round teeth
- W14 Partially erupted

MOLAR – WINGLESS (CONTINUED)

- W56 Large molar
- 26N Universal upper
- 28N General purpose lower molar

LABIAL – WINGED

- 6 Central and incisor
- 9 Labial prep, all anteriors

Tiger Clamps

- | | |
|------------------------------------|-------------------------------------|
| 1T Premolar, general purpose upper | 9T Anterior labial |
| 2T Premolar, general purpose lower | 14T Molar, partially erupted |
| 2AT Premolar, large neck premolar | 14DT Molar, distal access, serrated |
| | 56T For large molars |

Accessories

- | | |
|---|------------------------------------|
| Ivory Rubber Dams – green, blue, black and purple, 5" x 5" or 6" x 6" | Ivory Rubber Dam Punch |
| Ivory Rubber Dam Template | Ivory Rubber Dam Clamp Starter Kit |
| Ivory Rubber Dam Frames | Ivory Rubber Dam Starter Bundle |
| Ivory Rubber Dam Forceps | |